# Withings Thermo™

## Smart temporal thermometer

### Quick Installation Guide

---

**FR** Guide d'installation    **ES** Guía de instalación  
**DE** Installationsanleitung    **IT** Guida all'installazione

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# How to take a temperature measurement

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
**IT** Come misurare la temperatura

# 1

## Press the button to turn the thermometer on.

- FR Appuyez sur le bouton pour allumer le thermomètre.
- DE Drücken Sie auf die Taste, um das Thermometer einzuschalten.
- ES Presione el botón para encender el termómetro.
- IT Premere il pulsante per accendere il termometro.
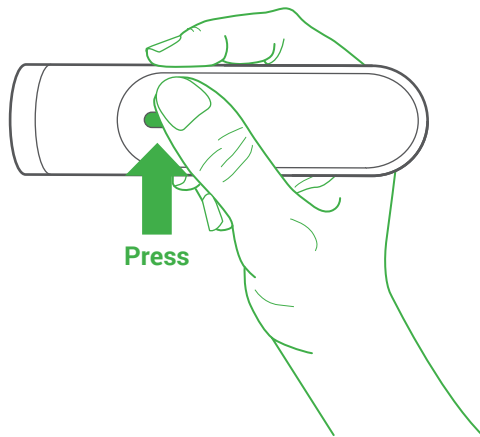
## 2

Starting from the middle of the forehead, scan straight across to the hair.

- FR En partant du milieu du front, balayez en ligne droite jusqu'aux cheveux.
- DE Ausgehend von der Mitte der Stirn, parallel zum Haaransatz messen.
- ES Explora a través de la frente, directamente desde el centro de la frente hasta el cabello.
- IT Passare sulla fronte, a partire dal centro della fronte fino all'attaccatura dei capelli.

### Clear any sweat or hair from the area before measurement.

- FR La zone de contact doit être dégagée et propre, sans sueur ni cheveux.
- DE Der Kontaktbereich muss frei und sauber sein, d. h., er darf weder von Schweiß noch von Haaren bedeckt sein.
- ES La zona de contacto debe estar limpia, despejada y sin sudor o cabello.
- IT La zona di contatto deve essere libera e pulita, senza sudore o capelli.

### 3


**The thermometer will vibrate when the measurement is complete.**

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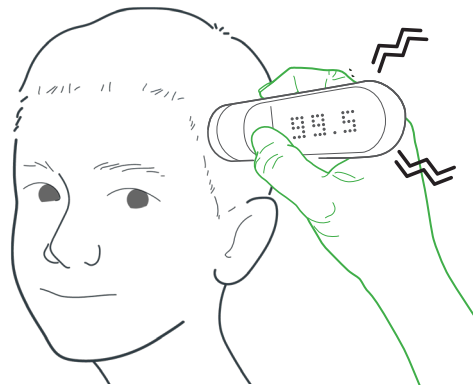# 4


**To assign the result, slide up and down on the end of Thermo to select a user (1). Press button to confirm your choice (2).**

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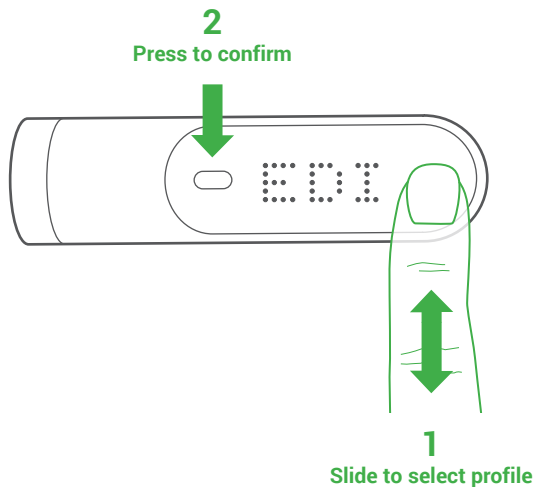# 5


## A colored LED indicates the fever level according to the user's age.

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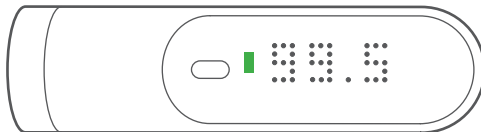
**High fever** | Fièvre élevée | Fieber | Fiebre alta | Febbre alta
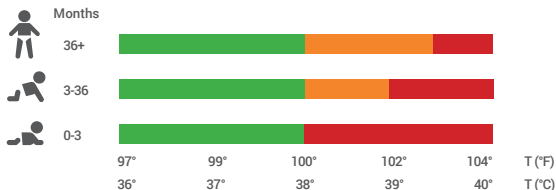
### Readings displayed by Withings Thermo are equivalent to rectal temperature levels

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- **Person being measured must be indoor for 15 minutes before taking a reading.**
- **Person being measured must have dry skin. Dry with a towel if needed before taking measurement.**
- **Person being measured and the thermometer should be in the same ambient temperature for at least 15 minutes before taking a reading.**
- **Person being measured should not have exercised over the last 15 minutes before taking a reading.**
- **Don't take a measurement while or immediately after nursing a baby.**

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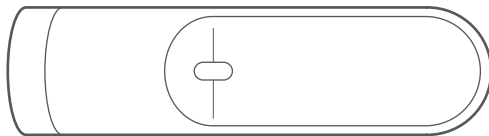
# Warranty

Garantie | Garantie | Garantía | Garanzia

## US WITHINGS One (1) Year Limited Warranty - Withings Thermo

Withings Inc. of 16192 Coastal Highway, Lewes, DE 19958 ("Withings") warrants the Withings-branded hardware product ("Withings Product") against defects in materials and workmanship when used normally in accordance with Withings' published guidelines for a period of ONE (1) YEAR from the date of original retail purchase by the end-user purchaser ("Warranty Period"). Withings' published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Withings does not warrant that the operation of the Withings Product will be uninterrupted or error-free. Withings is not responsible for damage arising from failure to follow instructions relating to the Withings Product's use.

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## UK WITHINGS One (1) Year commercial guarantee - Withings Thermo

Withings SA, 2 rue Maurice Hartmann 92130 Issy les Moulineaux, France ("Withings") warrants the Withings-branded hardware product ("Withings Product") against defects in materials and workmanship when used normally in accordance with Withings' published guidelines for a period of ONE (1) YEAR from the date of original retail purchase by the end-user purchaser ("Warranty Period"). Withings' published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Withings does not warrant that the operation of the Withings Product will be uninterrupted or error-free. Withings is not responsible for damage arising from failure to follow instructions relating to the Withings Product's use.

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**FCC ID: XNASCT01 IC ID: 11411A-SCT01**

This device complies with Part 15 of the FCC Rules and with Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation. This device complies with the R&TTE Directive. A copy of the EU Declaration of Conformity is available online at [Withings.com/compliance](http://Withings.com/compliance). **WARNING:** This product may contain one or more chemicals known to the State of California to cause cancer, birth defects, or other reproductive harms.

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©2016 Withings. All rights reserved.